SonoSite® Mini-Dock User Guide

The SonoSite® Mini-Dock provides connections for accessories, peripherals, and the M Series ultrasound systems, MicroMaxx® ultrasound system, or TITAN® high-resolution ultrasound system.

Attaching and removing the mini-dock

WARNING: To avoid the risk of electrical shock to the patient, do not touch the ungrounded connectors on the back of the mini-dock while simultaneously touching a patient.


Caution: To avoid damaging the mini-dock, do not perform high current grounding impedance testing involving this part. Accessible metal of the mini-dock is not protectively earthed. For an explanation on the grounding of the mini-dock, see SonoSite Technical Bulletin 003/06 at www.sonosite.com. Select **Support and Service** and then select your system type.

To avoid damaging equipment, remove mini-dock from system prior to transportation.

-
- Attach mini-dock**
- 1 Ensure that the ultrasound system is turned off.
Connecting the system and mini-dock while the power is on may cause the mini-dock to function improperly.
 - 2 Slide the mini-dock feet under the ultrasound system, and firmly press the mini-dock into the system connector. The mini-dock is correctly attached when no gap is between the system and mini-dock.

Note: If at any time the ultrasound system loses connectivity, turn off the system and reattach the mini-dock and system. Always verify that the mini-dock and system are correctly attached after changing a transducer. When you change a transducer, the mini-dock may loosen.

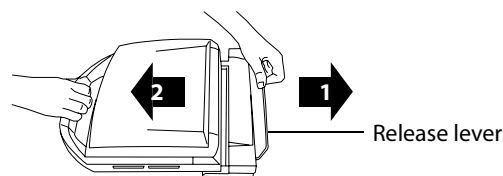-
- Remove mini-dock**
- 1 Grasp the handle on ultrasound system and push back on the release lever.
 - 2 Pull the ultrasound system forward.

Mini-dock connectors

Each connector on the mini-dock has a symbol that describes its use.

Connectivity Symbols on Mini-Dock

Symbol	Definition	Symbol	Definition
	DC input		S-video in
	Print control		DVI video out
	USB		Composite video out
	Ethernet		Audio out
	RS-232 (DVD recorder, VCR, or bar code scanner)		ECG/Footswitch
	S-video out		

Cleaning and Disinfecting

Caution: Use only recommended cleaners or disinfectants on surfaces. Immersion-type disinfectants are not approved for use on surfaces.

Clean and disinfect	<p>The exterior surface of the mini-dock can be cleaned and disinfected with a recommended cleaner or disinfectant.</p> <ol style="list-style-type: none">1 Remove the system, unplug the power supply, and detach any cables.2 Clean the surfaces using a soft cloth lightly dampened in a mild soap or detergent cleaning solution. Apply the solution to the cloth rather than the surface.3 Wipe the surfaces with the disinfecting solution. Theracide disinfectant is recommended.4 Air dry or towel dry with a clean cloth.
----------------------------	---

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Copyright 2007 by
All rights reserved.
SonoSite, Inc.
SonoSite Ltd